beurer

BY 90

DE	Babywaage Gebrauchsanweisung.....	2
EN	Baby scale Instructions for use.....	13
FR	Pèse-bébé Mode d'emploi	23
ES	Báscula para bebés Instrucciones de uso	34

IT	Bilancia pesaneonati Istruzioni per l'uso	46
TR	Bebek terazisi Kullanım kılavuzu	57
RU	Весы для младенцев Инструкция по применению	67
PL	Waga niemowlęca Instrukcja obsługi	78

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89231 Neu-Ulm

GermanyGermany

Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Dear customer,

Thank you for choosing a product from our range. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, gentle therapy, blood pressure/diagnosis, weight, massage and air.

With kind regards,
Your Beurer team

Table of contents

1. Getting to know your device.....	14	7. Additional information.....	20
2. Signs and symbols	14	8. Cleaning and maintenance	21
3. Intended use	15	9. Storage	21
4. Notes	15	10. Disposal.....	21
5. Device description.....	16	11. What if there are problems?	22
6. Initial use	17	12. Warranty/service	23

Included in delivery

- Brief instructions
- BY 90 baby scale
- 3 x 1.5V AAA batteries
- These instructions for use

1. Getting to know your device

A baby's physical development requires a great deal of energy. This is why regular weight control is important to put parents' minds at ease that their baby is getting sufficient nutrition and is developing healthily. With the specially developed baby scales from beurer, you will always be able to keep an eye on your child's weight to ensure their healthy development.

The baby scale has the following features:

- Curved weighing surface
- Large, easy-to-read LCD display
- Automatic and manual hold function for recording the weight of wriggling babies
- Integrated tape measure (150 cm)
- Tare weighing
- 20 kg weight capacity
- 5 g graduation
- Automatic switch-off
- Adjustable between weight units (kg/lb/oz)
- Storage of the last 10 measurements
- Transfer of measured values to the "beurer BabyCare" app
- The product uses *Bluetooth*[®] low energy technology, frequency band = 2400 MHz - 2480 MHz
Max. 5.0 dBm transmission power
Compatible with *Bluetooth*[®] 5.0 smartphones/tablets

System requirements

iOS from version 9.0, Android™ from version 5.0 with *Bluetooth*[®] ≥ 4.0.

List of compatible devices:

2. Signs and symbols

	WARNING Warning notice indicating a risk of injury or damage to health
	IMPORTANT Safety note indicating possible damage to the device/accessory
	Note Note on important information
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Dispose of packaging in an environmentally friendly manner
	Manufacturer
	This CE symbol certifies that the product complies with the requirements imposed upon the manufacturer by the European Community.

Certification symbol for products that are exported to the Russian Federation and members of the CIS

3. Intended use

The baby scale is only intended for measuring the weight of babies and small children. The baby scale is not a medical device. The baby scale is only intended for private use. It is not suitable for commercial or clinical use. The baby scale must only be used for the purpose for which it is designed and in the manner specified in the instructions for use. Any form of improper use can be dangerous. The manufacturer is not liable for damage resulting from improper or incorrect use.

4. Notes

WARNING

Safety notes

- Keep packaging material away from children (risk of suffocation).
- Swallowing batteries can be extremely dangerous. Keep the batteries and scale out of the reach of small children. Should a battery be swallowed, seek medical assistance immediately.
- Protect the baby scale from knocks, moisture, dust, chemicals, drastic changes in temperature and nearby sources of heat (ovens, heaters).
- Repairs must only be carried out by Beurer Customer Services or authorised retailers.
- Do not press the button violently or with pointed objects.

IMPORTANT

- Do not place any objects on the scale when not in use.
- The precision of the scale may be impaired by strong electromagnetic fields (e.g. induction hobs, mobile phones).
- Remove any transport locks.

General notes

- The scale's maximum capacity is 20 kg (44 lb, 704 oz). The results for weight measurement are displayed in 5 g increments (0.01 lb, 0.18 oz).
- In cases of differing measurements (between the scale and app), only use the measured values displayed on the scale.
- When supplied to the customer, the scale is set to weigh in "kg".
- We hereby confirm that this product complies with the European RED Directive 2014/53/EU. The CE Declaration of Conformity for this product can be found under:
<https://www.beurer.com/web/de/andingpages/cedeclarationofconformity.php>


WARNING

Notes on handling batteries

- If your skin or eyes come into contact with battery fluid, rinse the affected area with water and seek medical assistance.


- **⚠ Choking hazard!** Small children may swallow and choke on batteries. Store batteries out of the reach of small children.
- Observe the plus (+) and minus (-) polarity signs.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Protect batteries from excessive heat.
- **⚠ Risk of explosion!** Do not throw batteries into a fire.
- Do not charge or short-circuit batteries.
- If the device is not to be used for a relatively long period, take the batteries out of the battery compartment.
- Use identical or equivalent battery types only.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries.
- Do not disassemble, open or crush the batteries.

5. Device description

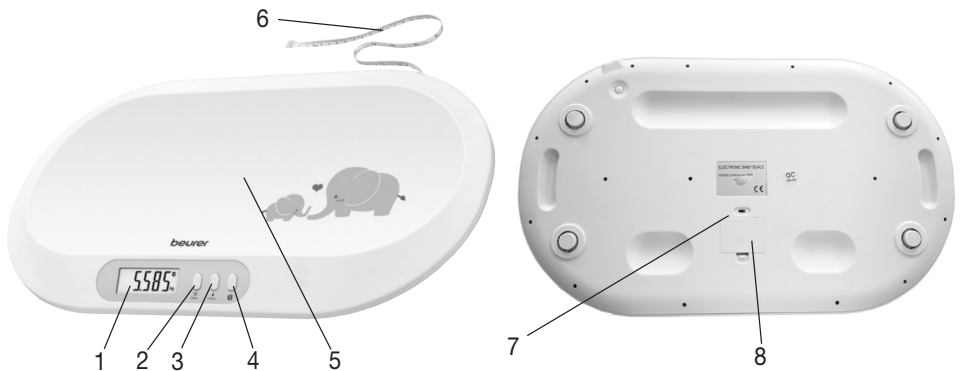1. Display	5. Weighing surface
2. ON/OFF/TARE button ①	6. Pull-out tape measure
3. HOLD button 🔒	7. Unit adjuster
4. Bluetooth® button 📶	8. Battery compartment with battery compartment cover

Display

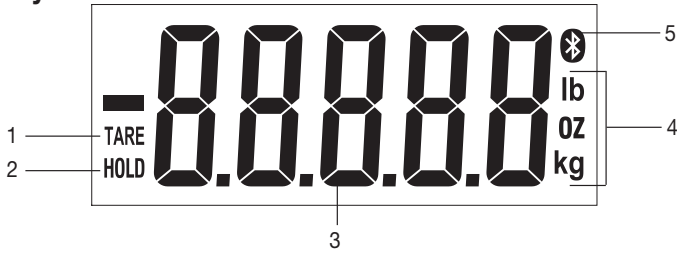1. "TARE" (TARE function active)	3. Weight display
2. "HOLD" (HOLD function active)	4. Weight unit (kg, oz or lb)
5. <i>Bluetooth</i> [®] symbol (<i>Bluetooth</i> [®] is activated)	

Selecting the unit of weight

When supplied to the customer, the baby scale is set to weigh in kilograms "kg". You can find the unit adjuster on the rear of the baby scale. You can use this to also change the unit of weight to pounds "lb" or ounces "oz".


6. Initial use

Inserting the batteries

Before using the baby scale, you must first insert three 1.5 V alkaline mignon AAA batteries into the battery compartment of the baby scale.

Proceed as follows:

1. Open the battery compartment cover on the bottom of the baby scale.
2. Insert three 1.5V alkaline mignon AAA batteries into the battery compartment. Make sure that the batteries are inserted the right way round (+/-), as shown in the battery compartment.
3. Close the battery compartment cover until you hear and feel it click into place.

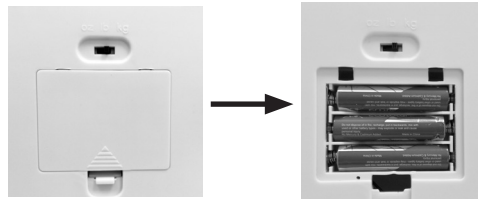i Note

The batteries must be replaced when "Lo" appears on the display.

Initial use with the app

You can transfer the measured weight values to your smartphone via *Bluetooth*[®] with the "beurer BabyCare" app.

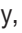
Keep the smartphone close to the scale in order to maintain an active *Bluetooth*[®] connection during initial use.

To connect the baby scale to the smartphone, proceed as follows:

1. Activate *Bluetooth*[®] in your smartphone's settings.

- ① Note: *Bluetooth*® devices, such as this scale, may not be visible in the general settings of the *Bluetooth*® device list.
- 2. Download the free "beurer BabyCare" app from the Apple App Store (iOS) or from Google Play (Android™).
- 3. Start the "beurer BabyCare" app and follow the instructions.
- 4. Enter user data.
The following settings must be set or entered in the "beurer BabyCare" App:

User data	Set values
Date of birth	DAY, MONTH, YEAR

- 5. Insert the batteries and place the baby scale on a solid surface ("app" will flash on the scale's display).
- 6. For 120 seconds the baby scale shows "app" in the display and can be connected to the "beurer BabyCare" app.
- ① Note: If you have not connected the baby scale with the app during initial use, "app" flashes for 10 seconds when the scale is switched on again. The indicator is no longer displayed once synchronisation with the app is successful.
- 7. The baby scale is successfully connected to the app once the *Bluetooth*® symbol  appears in the display.
- 8. The baby scale is successfully connected to the app once the *Bluetooth*® symbol illuminates continuously in the display.
- 9. Now place the baby on the weighing surface. Do this very carefully; do not let the baby fall and take particular care with the baby's head.
- 10. This baby scale has a 'HOLD' function (AUTO HOLD and MECHANICAL HOLD). The AUTO HOLD function (which is activated as standard) establishes the baby's weight on the display as soon as they are still enough to achieve an exact measured value. If the baby is moving around vigorously, their weight cannot be recorded. Therefore, try to calm the wriggling baby to achieve a rapid measurement, but do not touch the baby, as this will affect the weight measurement. As soon as the baby's weight is established on the display, the word 'HOLD' stops flashing and the weight is displayed for 120 seconds.
- ① Note: If the weight cannot be determined automatically after 8 seconds, the baby scale averages the weight over 4 seconds (flashing display of the weight on the baby scale display).
- 11. The weight is saved on the baby scale and it switches off automatically after 120 seconds. To switch off the scale manually, press and hold the ON/OFF/TARE button  for 3 seconds.
- 12. If the baby is wriggling about too vigorously, you can select the MECHANICAL HOLD function by pressing the HOLD button. After 3 seconds, the word 'HOLD' stops flashing and the averaged weight is displayed for 120 seconds.
- 13. To transfer the measured weight to the app, open the app and press the *Bluetooth*® button  after the measurement has successfully been taken.
- 14. The baby scale automatically switches off once the transfer has been completed successfully.

Initial use without the app

We recommend using the "beurer BabyCare" app for the initial use and implementing all settings (see "Initial use with the app" section).

i Note: Measurements that are saved on the baby scale without it having been previously synchronised with the app cannot be transferred.

1. Place the baby scale on an even, hard surface; a hard surface is crucial for correct measurement.
2. Press the ON/OFF/TARE button **I** to switch on the baby scale. "app" flashes on the display for 120 seconds. Press any button to use the baby scale without the app.
3. As soon as "0.000 kg" appears on the display, the baby scale is ready to take measurements.
4. Now place the baby on the weighing surface. Do this very carefully; do not let the baby fall and take particular care with the baby's head.
5. The weight is displayed for 120 seconds and then it switches off automatically. Press the ON/OFF/TARE button **I** to take another measurement. To switch off the scale manually, press and hold the ON/OFF/TARE button **I** for 3 seconds.

Subsequent synchronisation with the app

If you have not synchronised the baby scale with the "beurer BabyCare" app during initial use, you can also do this later.

1. Press the ON/OFF/TARE button **I** to switch on the baby scale.
2. "App" flashes on the display of the baby scale for 10 seconds. If the app is open, the baby scale can be connected for 10 seconds.

Alternatively:

Press the *Bluetooth*[®] button **B** whilst the baby scale is in operation. If the app is open, the baby scale can be connected for 120 seconds.

3. The baby scale is successfully connected to the app once the Bluetooth[®] symbol illuminates permanently in the display.

Transfer of the measurements

A maximum of 10 measurements can be stored.

- Press the *Bluetooth*[®] button **B**. **B** starts to flash in the display (120 seconds) until the app is opened. The *Bluetooth*[®] **B** symbol is then shown continuously in the display when there is a successful connection until the data is transferred. After successful transfer of the measurements, the scale switches off automatically.

Deactivating the AUTO HOLD function

If you wish to use the baby scale without the AUTO HOLD function, you can deactivate the AUTO HOLD function.

1. Once "0.000 kg" appears on the display, the AUTO HOLD function can be deactivated by pressing the HOLD button for 3 seconds.
2. As soon as "0.000 kg" appears on the display, the baby scale is ready to take measurements.
3. Now place the baby on the weighing surface. Do this very carefully; do not let the baby fall and take particular care with the baby's head.
4. To establish the baby's weight, press the HOLD button. The recording of the weight is confirmed in the display when the HOLD symbol flashes for 5 seconds.

5. The weight is displayed for 120 seconds and then it switches off automatically. Press the ON/OFF/TARE button **(I)** to take another measurement. In order to manually switch off the baby scale, press and hold the ON/OFF/TARE button **(I)** for 3 seconds.
6. To reactivate the AUTO HOLD function, press and hold the HOLD button for 3 seconds.

(i) Note: If the HOLD symbol + ON appear in the display, the AUTO HOLD function is activated.

TARE function

If you prefer not to place the baby directly onto the weighing surface, but instead onto a hand towel, for example, you must place the hand towel on the baby scale before switching it on. If you have placed a towel or other objects on the baby scale when it is already switched on and you would like to remove these again before weighing, you can also use the TARE function. Once you have pressed the TARE button, "TARE" is shown in the display. To reset the baby scale to "0.000 kg", you have to remove the towel or additional weight from the baby scale and press the TARE button again.

Using the tape measure

Pull out the tape measure on the top of the baby scale if you want to measure the baby. In order to reel the tape measure back into the baby scale, press the round button next to the tape measure on the bottom of the baby scale. You can manually enter the height and head circumference of your baby in the "beurer BabyCare" app.

7. Additional information

Saving the measured values on the baby scale and transferring them to the app

A maximum of 10 measurements can be saved in the baby scale. The saved measurements are transferred automatically to the app if you open the app within the *Bluetooth*[®] range and *Bluetooth*[®] is switched on on the scale.

It is not possible to automatically transfer data when the baby scale is switched off.

Deleting scale data

In order to completely delete all measurements on the baby scale, press and hold the *Bluetooth*[®] button **(B)** for 5 seconds.

Changing the batteries

Your scale is equipped with a low battery indicator. If you operate the scale with flat batteries, "Lb" will appear on the display and the scale will automatically switch off. In this case, the batteries must be replaced (3 x 1.5 V AAA). All saved measurements remain saved.

(i) Note:

- You have to re-synchronise the baby scale with the app (see "Initial use with the app" section) after replacing the batteries.
- When replacing the batteries, use batteries of the same type, make and capacity.
- Do not use rechargeable batteries.
- Use batteries free from heavy metals.

8. Cleaning and maintenance

The accuracy of the measurements and service life of the device depend on its careful handling:

The baby scale should be cleaned from time to time.

This should be done using a damp cloth and, if necessary, a small amount of detergent.

IMPORTANT

- Never use harsh solvents or cleaning products!
- Do not submerge the baby scale in water.
- Do not rinse the baby scale under running water.
- Do not clean the baby scale in the dish washer.

9. Storage

If you do not plan to use the baby scale for a long period of time, we recommend that you store it in the original packaging in a dry environment. During storage, do not place any objects on the baby scale. Store the baby scale away from children and pets. Remove the batteries from the baby scale.

10. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life.

Dispose of the device at a suitable local collection or recycling point in your country. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

Battery disposal

Empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

11. What if there are problems?

Display	Cause	Solution
Err	1. No flat, stable surface. 2. The maximum weight capacity of 20 kg has been exceeded.	1. Place the baby scale on a flat surface without carpet. Re-tare the baby scale. 2. The maximum weight capacity of 20 kg has been exceeded.
Incorrect weight is displayed.	No flat, stable surface.	Place the baby scale on a flat surface without carpet. Re-tare the baby scale.
Incorrect weight is displayed.	The baby is wriggling too much.	Try to calm the baby so that it stops wriggling.

Incorrect weight is displayed.	Incorrect zero point of scale.	Switch the baby scale off and on again. Wait until the baby scale displays "0.000 kg" and repeat the measurement.
No Bluetooth® connection (⊗ symbol is missing).	Device outside the range.	The range in the open is Approx. 15 m. Walls and ceilings reduce the range. Other radio waves may interfere with the transmission. Therefore, do not position the baby scale near devices such as WiFi routers, microwaves and induction hobs.
	Bluetooth® is switched off on the scale.	Press the <i>Bluetooth</i> ® button  until  lights up in the display.
	No connection to the App.	Completely close the App (including in the background). Switch <i>Bluetooth</i> ® off and back on again. Switch the smartphone off and back on again. Remove the batteries from the baby scale for a short period and reinsert them. See FAQs at www.beurer.com .
Lo	The batteries in the scale are empty.	Replace the batteries in the baby scale.
APP	The baby scale has not been synchronised with the App. All measurements carried out prior to synchronisation of the baby scale with the App cannot be transferred.	Connect the baby scale to the App (see "Initial use with the App" section).

12. Warranty/service

In case of a claim under the warranty please contact your local dealer or the local representation which is mentioned in the list "service international".

In case of returning the unit please add a copy of your receipt and a short report of the defect.

The following warranty terms shall Apply:

1. The warranty period for BEURER products is either 5 years or – if longer – the country specific warranty period from date of Purchase.
In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Repairs (complete unit or parts of the unit) do not extend the warranty period.
3. The warranty shall not be valid for damages because of
 - a. improper treatment, e.g. nonobservance of the user instructions.
 - b. repairs or tampering by the customer or unauthorised third parties.
 - c. transport from the manufacturer to the consumer or during transport to the service centre.
 - d. The warranty shall not be valid for accessories which are subject to normal wear and tear (cuff, batteries etc.).
4. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

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The *Bluetooth*[®] word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Beurer GmbH is under license. Other trademarks and trade names are those of their respective owners.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.

Android is a trademark of Google LLC.