
beurer

GS 39

- (D) Sprechende Glaswaage**
Gebrauchsanleitung
- (F) Pèse-personne vocal en verre**
Mode d'emploi
- (GB) Speaking glass scale**
Instruction for Use
- (E) Báscula de vidrio con voz**
Instrucciones para el uso
- (RUS) Стекланные весы**
Инструкция по применению

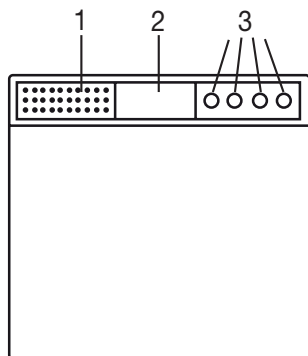(D) Gerätebeschreibung

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Dear Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for the applications in the areas of heat, weight, blood pressure, body temperature, pulse, gentle therapy, massage and air.

Please read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

Best regards,
Your Beurer Team

1. Interesting facts

Weighing correctly

If possible, always weigh yourself at the same time of day (ideally in the morning), after going to the toilet, on an empty stomach and without clothes to achieve comparable results.

2. Important notes – keep for later use!

■ Caution – For your safety

- Attention! Do not step onto the scale with wet feet or when the scale's surface is damp -danger of slipping!
- You should protect the scale against impacts, moisture, dust, chemicals, major temperature fluctuations and keep it away from heat sources (ovens, radiators).
- Do not stand on the outer edge to one side of the scale otherwise it may tip!
- Repairs may only be performed by Customer Service or authorized dealers. However, always check the batteries and replace them if necessary prior to making a complaint.
- All scales comply with EC Directive 2004/108 + supplements. Should you still have questions regarding the use of our products, please contact your dealer or Customer Service.
- The scale should occasionally be cleaned with a damp cloth. Do not use aggressive cleaning agents and never submerge the scale.
- The device is only for private use and is not intended for medical or commercial purposes.
- Please dispose of the appliance in accordance with the Directive 2002/96/EC on Waste Electrical and Electronic Equipment (WEEE). In case of queries, please contact the municipal authorities responsible for waste disposal in your area.

3. Voice function

The device has 5 language settings:

German (DE), French (FR), English (GB), Spanish (ES) and Russian (RUS)

The device is set to German when supplied to the customer.

To change the language, switch the scale on as described in the "4.2 Measuring weight" chapter and wait until "0.0" is shown on the display. Now press the middle button, "Language", on the back of the scale, until the desired language is announced.

Volume control:

You can adjust the volume of the language function as you wish, or switch it off by repeatedly pressing the "🔊 +/-" button on the back of the scale.

Note: Language off (no beep) means that voice output is not selected and the measurement value will not be announced.

4. Use

This scale is used to determine the weight of human beings. The scale has a weighing range up to 330 lbs (150 kg/24st).

4.1 Batteries

If present, pull the battery insulating strip off the battery compartment cover or remove the battery's protective film and insert the battery according to the polarity. If the scale fails to operate, remove the battery completely and insert it again. Your scale is equipped with a "replace battery" indicator. When the scale is operated with a battery which is too weak, "Lo" appears in the display field and the scale automatically switches off.

Display	DE	FR	GB	ES	RUS
Lo	Batterie schwach	Batterie faible	Battery weak	Pila débil	Слабый заряд батареек

The battery must be replaced in this case (3 x AAA 1,5V). The used, completely drained standard and rechargeable batteries must be disposed of in specially marked collection containers, at toxic waste collection points or electrical product retailers. You are legally obligated to dispose of the batteries.

Note: Batteries containing pollutants are marked with the following symbols: Pb = Battery contains lead, Cd = Battery contains cadmium, Hg = Battery contains mercury.

Place the scale on a firm, flat surface; a firm floor covering is a basic requirement for correct weight measurement.

When supplied to the customer, the scales are set to weigh in "kg". To the rear of the scales, there is a button „kg/lb/st" which you can use to change to 'pounds' and 'stones' **when the scale is in English.**

4.2 Measuring weight

Place the scale on a firm level floor (no carpet); a firm floor covering is required for correct measurement.

Audibly, briefly and firmly tap the scale with your foot to activate the vibration sensor!

Fig. 1

The entire display (Fig. 1) appears up to "0.0" (Fig. 2) as a self-test (tap-on technology).

Fig. 2

Display	DE	FR	GB	ES	RUS
8888	Beep	Beep	Beep	Beep	Звуковой сигнал

Now the scale is ready to measure your weight.

Display	DE	FR	GB	ES	RUS
0.0	Das Gerät ist betriebsbereit	L'appareil est prêt à fonctionner	The device is ready for use	El aparato está preparado para el servicio.	Прибор готов к работе

Now step onto the scale. Stand still on the scale with your weight distributed evenly between both legs. The scale immediately begins to measure your weight. After flashing twice the final measurement result appears on the display. (Fig. 3)

There is also an acoustic output. The scale switches off 10 seconds after you step down from it. The measuring result remains visible during this time.

Fig. 3

Fig. 4

If the scale detects an error during weighing, "----" or "Err" appears in the display (Fig. 4+5).

Fig. 5

Display	DE	FR	GB	ES	RUS
---- / Err	Fehler	Erreur	Error	Error	Неисправность

If a weight of more than 150 kg is placed on the scale, "OVER" is displayed (Fig. 6).

Fig. 6

Display	DE	FR	GB	ES	RUS
0/Err	Überlastung	Surcharge	Overload	Sobrecarga	Перегрузка

Make sure that you always switch on the scale first and wait for the display "0.0" (Fig. 2) before stepping on.

4.3 Weighing with storage function

The scale has 4 storage locations (buttons 1-4).

Switch on the scale as described in "4.2. Measuring weight" and wait until the display shows "0.0" (Fig. 2).

Now press the desired storage location (1, 2, 3 or 4).

For example, press the "1" button that appears in the upper left-hand corner of the display "User 1".

Display	DE	FR	GB	ES	RUS
User 1	Speicherplatz 1	Emplacement de mémoire 1	Memory space 1	Espacio de memoria 1	Пользовательски элемент позиция 1

Now the scale is ready to measure your weight, as described in 4.2.

The measurement is automatically saved in the respective user memory. If necessary, a previous measurement is exceeded.

4.4 Retrieving the stored values

Press the button for the desired storage location when the scale is switched off. The saved measurement is shown.

Display	DE	FR	GB	ES	RUS
User 1 62,4 kg	Speicherplatz 1 Benutzerspeicher 62,4 kg	Emplacement de mémoire 1 Mémoire 62,4 kg	Memory space 1 user memory 62,4 kg	Espacio de memoria 1 Memoria el usuario 62,4 kg	Пользовательски элемент позиция 1 позиция 62,4 kg

The complete display appears briefly as a self test (Fig.1), until "0.0" is displayed.

Display	DE	FR	GB	ES	RUS
0.0	Das Gerät ist betriebsbereit	L'appareil est prêt à fonctionner	The device is ready for use	El aparato está prepara- do para el servicio.	Прибор готов к работе

Now the scale is ready to measure your weight, as described in 4.2.

The new measurement is automatically saved in the selected storage location.

Tip: An other storage location can be only selected when "0.0" appears on the display.

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