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- If the device is not used for a long period, take the batteries out of the battery compartment.
- Use identical or equivalent battery types only.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries.
- Do not disassemble, split or crush the batteries.

#### 2. Preparation

In order to use the Callus remover, you first need to insert the supplied batteries. To insert the batteries, follow the steps below:

- Remove the battery compartment lid on the rear of the device.
- Insert the two supplied 1.5 V AA (alkaline type LR6) batteries. Make sure that the batteries are inserted with the correct polarity, according to the label. Do not use rechargeable batteries.
- Close the battery compartment lid carefully.


#### 3. Application

##### 3.1 General

- This device is intended for the treatment of feet only (pedicure).
- To change an attachment, press the unlock button and take the attachment off the device.
- If you wish to remove highly callused areas, use the coarse callus attachment.
- If you wish to remove slightly callused areas or to smooth your skin, use the fine callus attachment.
- Switch on the device by moving the sliding switch forwards. When doing so, also press the switch-on button on the sliding switch downwards.
- The sandpaper attachments are characterised by their extremely high effectiveness.
- Do not exert heavy pressure and always guide the attachments carefully over the surface to be treated.
- Move the instrument under light pressure in circular movements slowly over the parts to be treated.
- Note that the cutting attachments will have a reduced effect if the skin is softened or moist.
- For this reason, do not carry out any pretreatment in a bath.
- Do not remove all the callus skin in order to retain the natural protection of the skin.
- After every use, apply moisturising cream to the parts treated.

##### Changing the roller

- Proceed as follows to change the roller:
- Hold the callus remover in your hand as shown.
- With your other hand, push the roller to the left and pull it up and off.

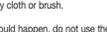- Clean the instrument with a slightly moistened cloth. You can also moisten the cloth with mild soap suds to remove more stubborn dirt.
- The sandpaper attachments must only be cleaned using a dry cloth or brush.
- Do not drop the instrument.
- Make sure that no water enters the instrument. But if this should happen, do not use the instrument again until it is completely dry.
- Do not use chemical detergent or scouring agents for cleaning.
- Never immerse the device in water or other liquids.

#### 5. Replacing the battery

If the inserted batteries are empty, you must replace them. To replace the batteries, follow the instructions in the "Initial use" chapter.

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#### 6. Replacement parts and wearing parts

Replacement part	Item number and/or order number
Sandpaper attachment replacement set	573.07

Commercially available or from the corresponding service address (according to the service address list)

#### 7. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the unit at a suitable local collection or recycling point. Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

The empty, completely flat batteries must be disposed of through specially designated collection boxes, recycling points or electronics retailers. You are legally required to dispose of the batteries.

The codes below are printed on batteries containing harmful substances:

Pb = Battery contains lead

Cd = Battery contains cadmium,

Hg = Battery contains mercury.